



WHO WE ARE

Pharma Machinery & Spares Company, located in Ahmadabad, Gujarat manufacturer & supplier of Pharma Machinery & Change Parts of Compression Machines, Tableting Equipment's Allied Machines of internationally accepted standards.

Our mission is to provide turnkey solutions to quality conscious global clients. We have approved list of vendors who will manufacture Pharma Machinery & Change parts of Compression Machine and Tableting Equipment's / Allied Machines as per our standards & specifications. We help our vendors to improve technology and product features.

We pride ourselves on consistently applying due diligence, high standards & excellent understanding of the value of methodology.

We are committed to quality maintenance, timely delivery schedules and adhering to buyers specifications.

We are exporting **Pharma Machinery & Change parts of** Compression Machine and Tableting Equipment's / Allied Machines to the most demanding market of Asian Countries.

The list of products, we have offered here is only indicative and we are always adding more products.

Range of Tablet Presses

- Single Rotary Mini Tablet Press / Compression Machine Rotary
- Single Rotary Tablet Press / Compression Machine with 16D/20D/23DB / 26 BB Stations
- Double Rotary Tablet Press / Compression Machine 27D / 27B / 35B / 45BB Stations
- Double Rotary High Speed Tablet Press/Compression Machine with 37D / 45B / 55BB / 61BB Stations
- Double Rotary High Speed Tablet Press/Compression Machine with Pre-compression
- Single Rotary High Speed Tablet Press/Compression Machine with Pre-compression
- Central Drive Rotary Tablet Press / Compression Machine (Pre-compression or Double Rotary Options)
- Slugging Tablet Press / Bolus Compression Machine for Veterinary

Tableting Equipment's/Allied Machines

- Tableting Tooling's Punch Die Sets for Tablet Press
- Dust Extractor Machine
- De Burring De Dusting Machine
- Oscillating Granulator
- Multi Mill
- Colloid Mill
- Coating Machine
- · Communiting Mill
- Sieving Grading Straining Machine

_

ROLL COMPACTOR MACHINE



LAB MODEL - Roll Compactor cGMP (Plain & Water Jacketed Model

| Model | CPM-LAB RC - 25 |
|--------------|-------------------|
| Roll Speed | 2-29 RPM Variable |
| Output * | 2 to 5 Kgs./Hr. |
| Screw Speed | 5-86 RPM Variable |
| Roller Motor | 3 HP / 3 Ph. |
| Feeder Motor | 2 HP / 3 Ph. |



ROLL COMPACTOR cGMP (Commercial)

| Model | CPM-LAB RC - 25 | CPM-LAB RC - 75 | CPM-LAB RC -200/100 | |
|------------------------|--------------------|-------------------------|----------------------|--|
| Output | 20 to 60 Kgs./Hr.* | 75 to 100 Kgs./Hr.* | 100 to 250 Kgs./Hr.* | |
| Roll size | Ø 200 mm x 50 mm | Ø 200 mm x 75 mm | Ø 200 mm x 100 mm | |
| Roll speed | 5 to 25 RPM | 5 to 25 RPM | 5 to 25 RPM | |
| Feed screw speed | 10 to 60 RPM | 10 to 60 RPM | 10 to 60 RPM | |
| Roll drive Motor | 5Hp/960RPM/3 Ph. | 5 Hp/960RPM/3 Ph. | 7.5Hp/960RPM/3 Ph. | |
| Feed screw drive motor | | 3 HP / 1440 RPM / 3 Ph. | | |



ROLL COMPACTOR cGMP

| Model | CPM-RC-200/50 | CPM-RC- 200/75 | CPM-RC- 200/100 |
|------------------------|--------------------|-------------------------|----------------------|
| Output | 20 to 60 Kgs./Hr.* | 75 to 100 Kgs./Hr.* | 100 to 250 Kgs./Hr.* |
| Roll size | Ø 200 mm x 50 mm | Ø 200 mm x 75 mm | Ø 200 mm x 100 mm |
| Roll speed | 5 to 25 RPM | 5 to 25 RPM | 5 to 25 RPM |
| Feed screw speed | 10 to 60 RPM | 10 to 60 RPM | 10 to 60 RPM |
| Roll drive Motor | 5Hp/960RPM/3 Ph. | 5 Hp/960RPM/3 Ph. | 7.5Hp/960RPM/3 Ph. |
| Feed screw drive motor | | 3 HP / 1440 RPM / 3 Ph. | |

- ➤ Totally cGMP Compact design availably in plain water jeckted & Flame proof model.
- ➤It's improve the bulk density of material.
- ➤It's improves uniform mixing of material.
- ➤ It's achieve better granulation of sieves analysis.
- ➤ It's produce uniform particle size range of granules.
- ➤ Various types of roll available to achieve maximum bulk density of out put.
- ➤ Specially designed feeding system consisting of hopper pre-densifer & variable speed drive unit.
- > Feeding unit can hydraulically lift by hydraulic pump for easy cleaning.
- ➤ Specially design provide in the roller to prevent powder sticking on the roll.
- ➤ Pre-compression chamber guides the powder without any slippage or leakage.
- >To prevent major accident the roll is designed to rotate reverse & forward when cleaning. Specially limit switch is provided to prevent major accident.
- Main roll drive is design with twin shaft gear box coupled with coupling to protect roll drive unit there fore perfectly synchronizing the roll with a constant torque

GRANULATION & COUNTING MACHINE



DOUBLE CONE BLENDER cGMP

| Model | Gross Volume (Ltr.) | Working Volume (Ltr.) | Motor (HP) |
|-------------|---------------------|-----------------------|------------|
| CPM-DCB-5 | 15 | 9 | 0.5 |
| CPM-DCB-10 | 30 | 18 | 0.5 |
| CPM-DCB-15 | 45 | 27 | 0.5 |
| CPM-DCB-25 | 75 | 45 | 1 |
| CPM-DCB-50 | 150 | 90 | 2 |
| CPM-DCB-100 | 300 | 180 | 3 |
| CPM-DCB-150 | 450 | 270 | 3 |
| CPM-DCB-200 | 600 | 360 | 5 |
| CPM-DCB-250 | 750 | 450 | 5 |
| CPM-DCB-300 | 900 | 540 | 7.5 |
| CPM-DCB-350 | 1050 | 630 | 7.5 |
| CPM-DCB-400 | 1200 | 720 | 10 |
| CPM-DCB-500 | 1500 | 900 | 10 |



RIBBON BLENDER cGMP

| _ | | _ | |
|------------|--------------------|------------|-----|
| Model | Total Cap(Ltr.) | R.P. M. | H.P |
| CPM-RB-20 | 30 | 100 | 0.5 |
| CPM-RB-60 | 90 | 67 | 1.0 |
| CPM-RB-200 | 300 | 46 | 3.0 |
| CPM-RB-300 | 430 | 40 | 5.0 |
| CPM-RB-500 | 700 | 34 | 7.5 |
| | | | |



TABLET COUNTER MACHINE cGMP

| Model | CPM - TC |
|----------------------|-----------------------|
| Tablet size | 5mm to 20mm |
| Motor | 0.5 HP/1440 RPM/3 Ph |
| RPM | 5 to 10 RPM |
| Net Weight(Aproxx) | 300 Kgs |
| Gross Weight(Aproxx) | 550 Kgs |
| Machine Dimension | 915L x 990B x 2515H |
| Grose Dimension | 1065L x 1140B x 2665H |



RAPID MIXER GRANULATOR cGMP

SALIENT FEATURES :

- >Models available from 15 Ltr to 800 Ltr working capacity
- Fully automatic operation based on PLC control
- ➤ Better productivity & lesser mixing time
- Contact parts are SS 316/316L
- ➤ Uniform mixing & granulation
- ➤ Less handling, more hygienic
- >cGMP with documentation

RAPID MIXER GRANULATOR cGMP



| Model | CPM - OGR | Model | CPM - OGR |
|--------------------|-----------------------------|--------------------|-----------------------------|
| Hopper capacity | 33 Liters | Hopper capacity | 33 Liters |
| Screen area | 388 mm (W) x 762 mm (L) | Screen area | 388 mm (W) x 762 mm (L) |
| Stirrer size | Ø 138 mm x 368 mm Length | Stirrer size | Ø 138 mm x 368 mm Length |
| Oscillation | 110 Per Minute. | Oscillatio n | 110 Per Minute. |
| Output * | 30 to 90 Kgs./Hr. | Output * | 30 to 90 Kgs./Hr. |
| Charging height | 1135 + 5 mm | Charging height | 1135 + 5 mm |
| Discharging height | , | | 655 + 5 mm |

MODEL AVAILABLE FROM SINGLE STIRRER & DOUBLE STIRRER

- ➤ Sturdy and noiseless design.
- ➤ Easy to operate & maintain.
- ➤ Higher ratio of desired size granules can be obtained
- ➤ Continuous oscillating stirrer drive.
- ➤ Simple screen holding and tensioning device.

TABLETING MACHINE

R&D LAB PRESS cGMP



| Model | CPM-LAB-8 D | CPM-LAB-10 B | CPM-LAB-12 D-B-BB | |
|-------------------------|------------------------------------|--|-------------------|--|
| Number of Stations | 8 | 10 | 12 | |
| Type of Tooling | D | В | D-B-BB | |
| Output - Tablets/Hr. | 2000 TO 14,000 | 2,000 TO 14,000 2,000 TO 18,000 2,000 TO | | |
| Turret RPM (Min. / Max) | 10 - 30. | 10 - 30. 10 - 30. | | |
| Max. Operating Pressure | 6 Tons. 6 Tons. 6 Tons. | | 6 Tons. | |
| Max. Tablet Dia. | 23 mm 16 mm 23 mm | | 23 mm | |
| Maximum Depth of fill | 17.5 mm 17.5 mm 17.5 mm | | | |
| Power | 2.0 HP X 1440 RPM X 415 V. X 50 Hz | | | |
| Overall Dimensions (mm) | 840 x 635 x 1200 H | | | |

- Imported needle roller bearing provided at the bottom of turret.
- >Upper punch penetration, as a standard features
- ➤ Electronic Digital Tablet cum RPM Counter, as a standard features.
- ➤ Guards of Acrylic Material.

- >Grease Nipples are provided for effective lubrication, where required
- ➤ Turret of Special Grade S.G. Iron casting
- ➤ Paint free tablet manufacturing zone
- ➤ Aluminum Turret Guards
- ➤ Square GMP Table Top Model
- ➤ Anti Vibating mount.



SINGLE ROATRY Pharma Machinery & Spares MACHINE cGMP

| Model | CPM-16 CPM-20 0 | | CPM-12 |
|-------------------------|--------------------------|-----------------------------------|--------------------------|
| Tooling | "D" | "D" E | |
| No. of Station | 16 Stn. | 16 Stn. 20 Stn. 1 | |
| Operating Pressure | 10 Ton | 10 Ton 10 Ton 1 | |
| Max. Tablet Diameter | 23 mm 23 mm 4 | | 40 mm |
| Max. Depth of Fill | 20 mm 20 mm 2 | | 25 mm |
| Output - Tablets/Hr. | 14,400 to 40,300 | 14,400 to 40,300 21,600 to 50,400 | |
| Upper Punch Penetration | 1.5 to 8 mm 1.5 to 8 mm | | 1.5 to 8 mm |
| Main Electric Motor | 2HP/960RPM/415V/50HZ/3Ph | | 5HP/960RPM/415V/50HZ/3Ph |

- ➤ Lower Guard of polished Stainless Steel & Upper Guards of Acrylic Material.
- ➤Inter lock switches to all Guards.
- ➤ Effective Dust Extraction Nozzles.
- ➤ Aluminum Turret Guards
- Reverse direction motion protection is provided to avoid accident.
- ➤ Imported needle roller bearing provided at the bottom of turret.
- ➤ Paint free tablet manufacturing zone.
- ➤ C.I. Body and C.I. Middle plate
- ➤ Turret of Special Grade S.G. Iron casting
- ➤ Upper punch penetration system.
- ➤ Electronic. Digital Tablet cum RPM Counter.
- ➤ Heavy duty square GMP Model.





DOUBLE SIDED ROTARY TABLETING MACHINE cGMP

| Model | CPM-D-27 | CPM-B-35 | CPM-BB-45 | |
|-------------------------|--|----------|--------------------|--|
| Type of Tooling | D | В | ВВ | |
| No. of Stations | 27 | 35 | 45 | |
| Max. Operating Pressure | 10 Tons | 6.5 Tons | 6.5 Tons | |
| Max. Tablet Diameter | 23 mm (25 mm on demand) | 16 mm | 11 mm | |
| Max. Depth of Fill | 20 mm 17.5 mm 17.5 mm | | 17.5 mm | |
| Output Tablets/Hr.* | 58,320 to 1,36,000 75,600 to 1,76,400 57,200 to 2,26,800 | | 57,200 to 2,26,800 | |
| Upper Punch Penetration | 1.5 to 8 mm 1.5 to 8 mm 1.5 to 8 mm | | | |
| Main Electric Motor | 5Hp / 1440 RPM / 415 V / 50 Hz / 3 Ph | | | |

➤ Heavy Rigid Body With Less Vibration.

➤ C.I. Body and C.I. Middle plate.

➤ Paint free tablet manufacturing zone.

➤ Square GMP Model.

➤ Turret of Special Grade S.G. Iron casting

> Aluminum Turret Guards.

➤ Double Sided lifting cams.

➤ Imported needle roller bearing provided at the bottom of turret.

➤ Effective Dust Extraction Nozzles.

➤ Variable Speed pulley, One shot Lubrication System, Anti Vibrating mount.

> Electronic. Digital Tablet cum RPM Counter. > Reverse direction motion protection is provided to avoid accident.



HIGH SPEED Pharma Machinery & Spares cGMP

| Model | CPM-Press-37D | CPM-Press-45B | CPM-Press-55BB | CPM-Press61BB |
|----------------------------|---------------------------|-----------------------|-------------------------|-------------------------|
| No. of stations | 37 | 45 | 55 | 61 |
| Type of tooling | D | В | BB | BB |
| Max. Operating Pressure | 10 Ton | 6.5 Ton | 6.5 Ton | 6.5 Ton |
| Max. Tablet Dia. | 23 mm (25.4 mmon demand) | 16 mm | 11 mm | 11 mm |
| Depth of Fill | 20.5 mm | 17.5 mm | 17.5 mm | 17.5 mm |
| Output Tablets / Hr.* | 79,920 to 2,22,000 | 97,200 to 2,80,000 | 1,18,800 to 3,83,000 | 1,31,760 to 3,80,640 |
| Upper Punch Penetration | 3.2 to 8 mm | | | |
| Main Electric Motor | 7.5 Hp / 1440 RPM / 3 Ph | | | |
| Powder Feeder Motor | 0.25 Hp / 1440 RPM / 3 Ph | | | |

- Front and easy outside controls.
- High Speed Pharma Machinery & Spares with outstanding engineering, performance and reliability.
- Auto lubrication with fine tuning valve and timer.
- AC invertor variable drive for main and feeder motor.
- ➤ Separate control panel safe against electric short Circuit during steam wash.
- Extremely robust quality engineering, house in an easily accessible cabinet.
- All the tablet parameters can be fine tuned during the operation of the machine from out side the tableting zone.
- Simple and convenient operation, operator protection, cleanliness and ease of maintenance.



ALLAID MACHINE



COLLOID MILL cGMP (Model Available from 15 Ltr to 35 Ltr)

| Model | CPM - CM-15 |
|---------------------------|----------------------------|
| Output * | 150 to 1500 Liters. |
| Particle size reduction * | 5 to 15 Microns. |
| Hopper capacity | 15 Liters |
| Main Electrical Motor | 3 HP / 2850 RPM / 3 Ph. |

- ➤ R&D Model Available on Customer Demand.
- ➤ Totally leak proof design.
- ➤ Unique for small volume micronization.
- ➤ Modified design of stator-rotor assembly.
- Available in plain, water jacketed & flame proof options.
- Easy to dismantle for cleaning and maintenance.



COATING MACHINE cGMP

| Model : CPM-CM | 12" | 18" | 30" | 36" | 42" | 48" |
|---------------------------|-----|-----|-----|-----|-----|-----|
| Loading Capacity in Ltrs. | 15 | 25 | 45 | 75 | 125 | 185 |
| Main Electric Motor | 0.5 | 0.5 | 2.0 | 2.0 | 3.0 | 3.0 |
| Blower Ph.0.5 | 0.5 | 1.0 | 1.0 | 1.0 | 1.0 | 2.0 |
| Heater Capacity in Kw. | 1.0 | 1.0 | 3.0 | 3.0 | 3.0 | 3.0 |

- ➤ Hot air blowing system is with temperature control probe.
- ➤ Coating Pan is engineered to the perfection from SS 304 (SS 316 optional).
- ➤ Totally enclosed machine drive and hot air blowing system.
- ▶12" or 48" Coating Machine has variable speed drive. Hot air blowing system is separate.
- >60" or 72" Coating Machine has two speeds with step pulley without hot air blowing system.




DE-BURRING & DE-DUSTING MACHINE cGMP

| Model | CPM - DDM |
|----------------------|-----------------------------------|
| Output * | 5,000 to 2,50,000 tablet / Hr. |
| Main Electric Motor | 0.25 Hp. / 1440 RPM / 3 Ph. |
| Max. Tablet Diameter | 25 mm |
| Min. Tablet Diameter | 5 mm |

- ➤ Components can be easily washed or sterilized.
- ➤ Designed for effective De-Burring and De-Dusting
- Noiseless & Maintenance free unit.
- Adjustable height to suit any tableting machine outlet.
- ➤ Minimum electric power consumption.
- Suction outlet concentrated on De-Dusting zone ensures maximum flushing of removed dust.

WE ARE ALSO IN

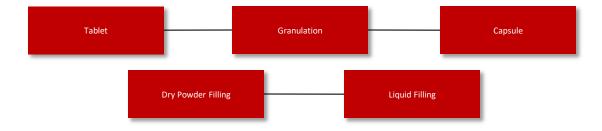



DUST EXTRACTOR MACHINE cGMP

| Model | CPM - DE-150 / 300 CFM |
|-----------------------|---|
| Electric Motor | 1 HP / 2800 RPM / 3 Ph. |
| Filtration Area | 1.343 Sq. Mtr. |
| Dust Storage Capacity | 40 Liters. |
| Suction Capacity | 150/300 CFM (500 Mtr. Cube / Hr.) at 125 mm W.G. |

Provided with stainless steel 304 contact parts (SS 316 optional).

Fabricated for effective dust collection.

Dynamically balanced blower fan ensures minimum noise

Easy to maintain and operate.

Four inlet manifolds facilitate to connect more dusting

Castor wheel provided for ease in mobility.

DRYING /PACKING MACHINE, SPARE PARTS



FLUID BED DRYER cGMP Available from 5 Kg. to 250 Kg.



TRAY DRYER cGMP Tray Dryer Available from 6 Tray / 12 Tray / 24 Tray / 48 Tray



TABLET INSPECTION MACHINE cGMP



UPPER CAM DISK





VARIABLE SPEED **PULLY**



FEED FRAME



TURRET



MULTI MILL SEIVE



SIFTER SEIVE



WORM



TRACKS



DIES & PUNCHES

CONTACT US



+91 7807506765 / +91 8320340951



sales.pharmamachinery@gmail.com

-HEAD OFFICE & FACTORY



Karnavti Industrial Estate, Nr. ONGC Well, Odhav, Ahmadabad – 38215 (Gujarat) India.

-BRANCH OFFICE & WORK SHOP



Suncity Road, Near Toll Plaza, Baddi, Dist – Solan, Himachal Pradesh - 173205 (INDIA)

WE ARE ALSO US ON

